



ENDOSCOPY - SURGERY - DENTAL



A collection of endoscopic instruments, including several long, thin shafts with various tips and handles, some with blue accents, arranged diagonally across the page. The background is a light gray gradient.

TEKNO LIZARD

TEKNO LIZARD

TEKNO LIZARD

TEKNO LIZARD

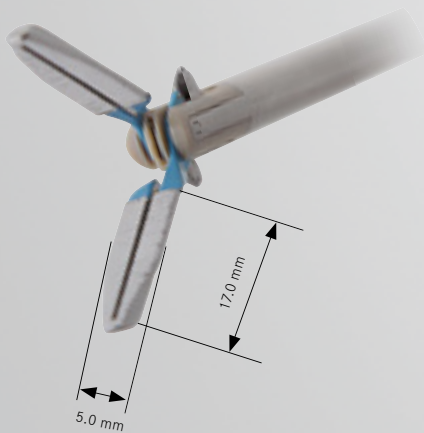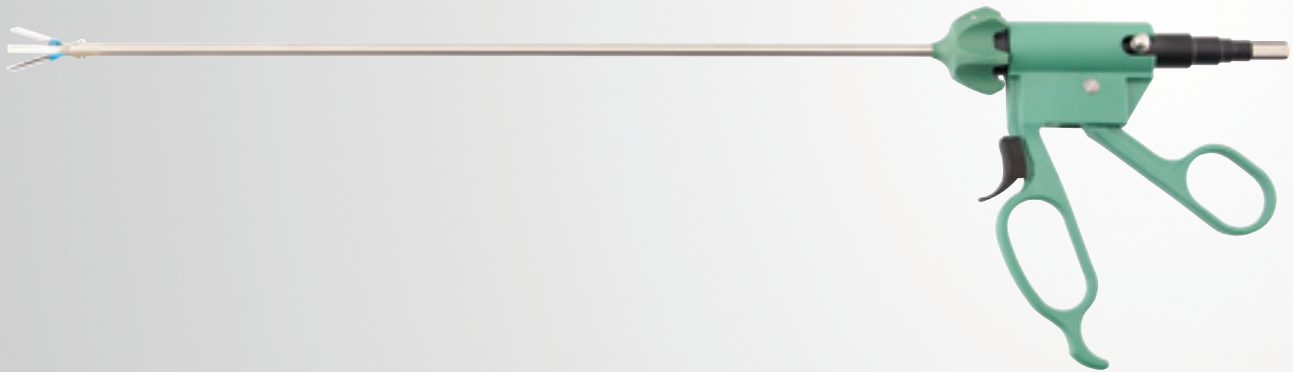

	250 mm	340 mm
Ø 5.0 mm	tk 707-060	tk 707-061

Tekno Lizard komplett, bestehend aus:
Tekno Lizard, complete, consisting of:

Ringgriff mit Standard-Flachstecker <i>Ring handle with standard flat plug</i>	1x	tk 707-050	tk 707-050
Außenschaft <i>Outer shaft</i>	1x	tk 707-051	tk 707-052
Einsatz (Elektrode) <i>Insert (Electrode)</i>	1x	tk 707-054	tk 707-053
Schubstange <i>Push rod</i>	1x	tk 707-057	tk 707-056
Halteschaft für Klingen <i>Blade holder shaft</i>	1x	tk 707-059	tk 707-058
Wechselbares Messer, unsteril, Pack á 20 Stück <i>Disposable blade, not sterile, box of 20 pcs.</i>	1x	tk 707-055	tk 707-055

Die Bauweise als Kombi-Instrument erlaubt das sichere Koagulieren von Blutgefäßen und Gewebe sowie anschließenden, mittigen Schnitt. Dies erspart dem Operateur den aufwändigen Wechsel der Instrumente sowie die Positionierung für den Schnitt. Die Wiederverwendbarkeit des Instruments erlaubt einen kostengünstigen Einsatz. Die Klinge wird nur einmal verwendet und lässt sich einfach und schnell einsetzen. Der Außenschaft ist 360° drehbar und ermöglicht somit dem Anwender die Elektrode nach seinen Wünschen zu positionieren.

The instrument combination allows the secure vessel coagulation followed by a center cut. Thus any change and re-positioning of instruments becomes obsolete. The instrument is reusable and cost efficient. The blade is quick and easy changed after each use. The outer sheath is 360° rotatable, permitting to place the electrode exactly according to the surgeon's needs.

© 2016 – TEKNO-MEDICAL Optik-Chirurgie GmbH
TEKNO is a registered trademark of TEKNO-MEDICAL Optik-Chirurgie GmbH
Printed in Germany 5/2016 – Nachdruck, auch auszugsweise, verboten